

# PREMIER Tags

## Instructions for Inserting Premier "X" Tags



### Fitting Instructions:

1. Place the female part of the tag completely under the spring fingers. The numbers will be down and the tag will be extending away from the tagger. Slide the male tag over the metal pin. The numbers will be up and the tag will extend away from the tagger.
2. Before you try to insert the tag into the ear, gently squeeze the applicator together to ensure that both the male and female tags are in alignment.
3. Place the ear between the jaws of the applicator, positioning the male tag on the back side of the ear. Make sure the tag is correctly placed as shown.
4. While gripping the ear quickly and firmly at the proposed position of the tag, place the point of the tag in contact with the ear close to your thumb. This allows you to feel the vein therefore avoiding it. Completely close the applicator until the tag is inserted through the ear and the female side of the tag.

### Important:


Use the black adaptor ONLY when inserting 2X or 5X Non-TamperProof Tags. Remove adaptor when inserting 2X or 7X TamperProof Tags.

### Tag placement is critical to success!

Tags, especially tags for sheep, must be installed no more than 2 inches from the skull and near the center of the ear.

The opposite extreme (installing too far from the skull) will result in more lost tags as the tissue on the tip of the ear is tender & tears easily.

Be careful to avoid the large vein (shown at right).

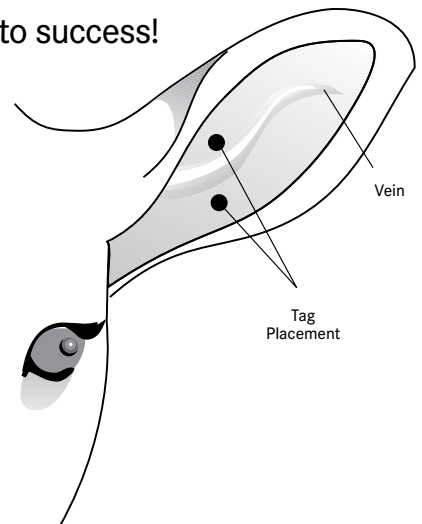

800-282-6631  
[www.premier1supplies.com](http://www.premier1supplies.com)  
Washington, IA